



Shakers

BIS Series - Incubated Shakers

BS Series - Orbital Shakers

being Introduction

BEING is an economically priced, high-end, high-performance laboratory equipment brand. We are committed to providing users with intelligent, intuitive, and professional laboratory equipment that modern laboratories require.

Besides the BS and BIS orbital shakers, BEING offers laboratories a broad portfolio of incubators, ovens, heated magnetic stirrers, water baths, circulators, chillers, and vacuum pumps.

BS / BIS Orbital Shakers



Our orbital shakers are the 'Smart Choice' for agitating and mixing samples and solutions and culture growth.

With 6 different models to choose from, BEING offers one of the largest selections of incubated and platform lab orbital shakers on the market. They are ideal for material dissolution, liquid mixing, cell culture, protein expression, fermentation, water quality analysis, DNA amplification, and many more applications in chemistry, clinical, pharmaceutical, and research laboratories.

All of our shakers are energy efficient, have excellent temperature regulation capabilities, and come with a host of features that provide safe and easy operation — and are economically priced. They're all designed, manufactured, and tested to the DIN 12880-2007 standard, providing a long service life.


This combination of selection, specifications, features, quality, and value makes BEING shakers the smart choice.



Precise temperature control


BEING BIS series incubated shaker provide a precise and steady heating environment that ensures consistent product quality, lowers the chances for rework and helps achieve reliable production results while reducing your laboratory's energy costs by being energy efficient.

Temperature Uniformity

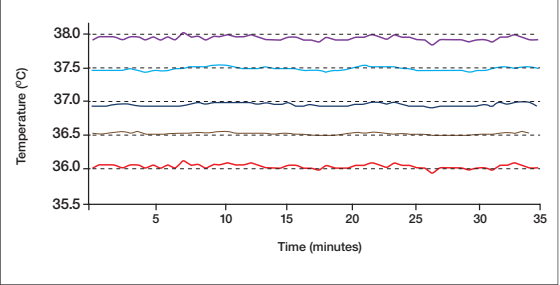

The incubated shaker's chamber's temperature uniformity ($\leq \pm 1.5^\circ\text{C}$) enables all samples to be heated evenly.


Temperature Stability

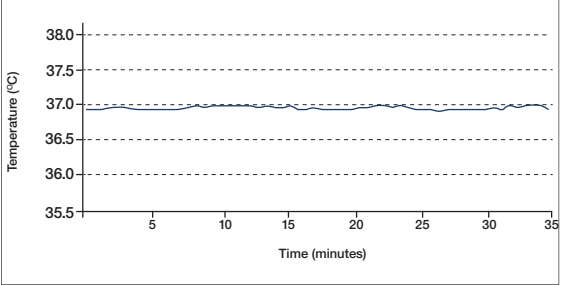

The incubation chamber's temperature stability of $\pm 0.1^\circ\text{C}$ ensures experiment stability.

Note: The stability and uniformity are measured at steady-state with an empty chamber according to DIN 12880.


Step Heating Control

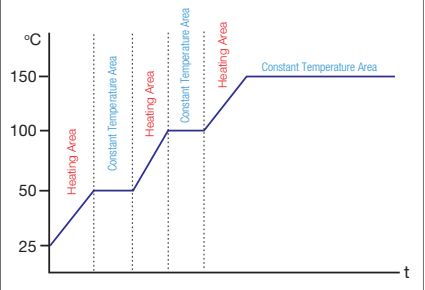


In program control mode, the controller allows the operator to set up step heating control.

Controller & Safety Feature		Incubated Orbital Shakers	Orbital Shakers
Series		BIS	BS
Controller	Automatic power on/off	√	√
	PID automatic control	√	√
	Programmable functions	√	√
	Fixed-value programs	√	√
	Multi-step programs	√	√
	Program cycling	√	√
	Timed & Untimed	√	√
	Oscillation speed - Adjustable	√	√
RUN delay	√	√	
Ports	Accessory Access	√	X
	Drain	√	X
Safety	Over-temperature protection	√	X
	Temperature limit protection	√	X
	Over-current protection	√	√
	Power off memory	√	√
	Anti-scalding protection	√	X
	Load Imbalance	√	√
	Audible & visual alarms	√	√

Professional LCD Controller Features

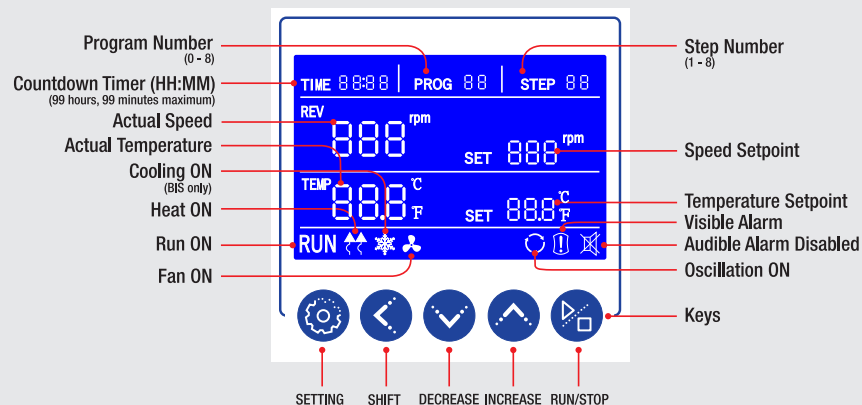


The BEING L-series professional controller has a 3.1-inch easy-to-read, real-time LCD display to select the various functions with all the parameters — temperature, time, program number, and step — on a single screen, providing quick and easy setting of temperature, time, and other parameters — and convenient operation.

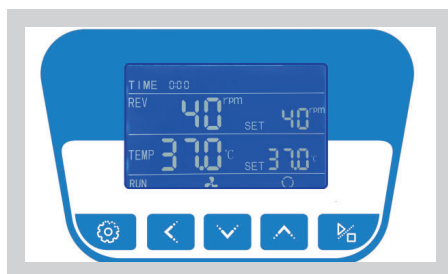


Easy to use

Simple and intuitive setting of all the operating parameters thanks to easy-to-understand icons and symbols.

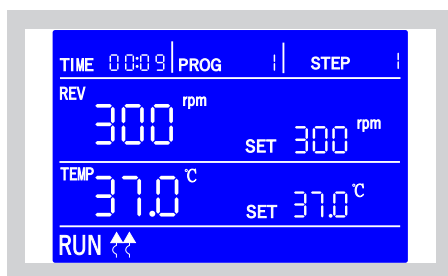

Programmable control at your fingertips


The shaker is designed to be used immediately out of the box with a single (fixed value) basic timed (1 minute to 99 hours and 59 minutes) or untimed program. But, if you need to simplify complicated testing processes and realize automatic control, that can be easily achieved. You can store and run up to 8 multistep programs with 8 steps each. Need to start your experiment at a specific time? The On/RUN delay makes it easy to do. Want to cycle your multistep program numerous times? That can be done up to 99 cycles.

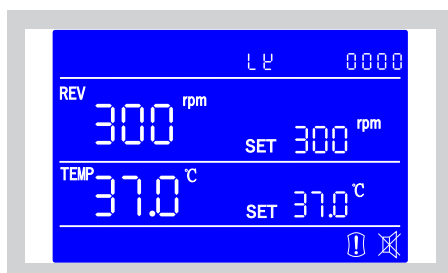

Password Protected

The controller has 4 settings access levels. Each level is password protected to avoid accidental changes to “sensitive” parameters.

Shaker Features

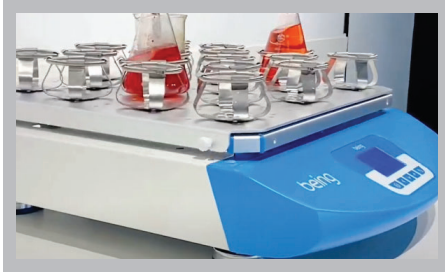

Toolless Accessory Platform Connection

The BEING orbital and incubated shakers are equipped with thumb screws that enable quick, easy, and toolless connection of platforms.



Patented Drive Technology

BEING shakers' rotation is controlled by a patented, maintenance-free, single-axis, eccentric drive and balancing technology utilizing a low energy consumption brushless DC motor that provides a large start torque and stable operation with a wide speed adjustment and low noise emission.



Large Transparent Chamber Door (BIS only)


The large transparent chamber door allows researchers and lab technicians to observe the incubation effects on the samples in real-time.



Gas Strut Door Support

Two gas struts smoothly control the opening of the BEING BIS Series incubated shaker's chamber door. They also hold the chamber door in any position. Platform rotation, heating element, and chamber fan are shut down when the chamber door is raised more than 25°.

Shaker Features

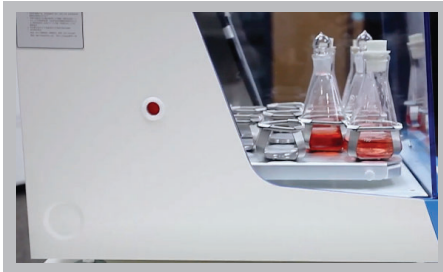

Temperature Test Hole (BIS only)

BEING BIS Series incubated shakers have an Ø18mm test hole and accessory port to allow multiple external temperature probes to be inserted into the chamber to validate temperature settings to the actual chamber temperature.



Condensate Drain (BIS only)

BEING BIS Series incubated shakers are designed to whisk away condensate that is generated during incubation. The chamber's bottom has channels on the edges to collect and direct the condensate to the drain on the back of the unit. The drain can be opened and closed with a push of a button.



Preventing damage from overheating (BIS only)

All incubators have dual overheating protection to prevent specimen and equipment damage. The controller's over-temperature protection is adjustable. It shuts down the heating element and fires an alarm until the temperature drops below the set point if the incubator exceeds the set temperature and the alarm setting. The independent overheating switch is adjustable.



Orbital Shaker


Model	BS-1	BS-2	BS-3
Maximum Capacity (ml)	2000	4000	9000
Platform Size (W x D) (in / mm)	10.0 x 10.0 250 x 250	13.8 x 13.8 350 x 350	17.7 x 17.7 450 x 450
Speed Range (rpm)	40 - 250	40 - 300	
Frequency Accuracy (rpm)	±1		
Amplitude (mm)	20		
Control	LCD Microprocessor PID		
Display Resolution	0.1		
Timer (hh:mm)	00:01 – 99:59		
Exterior Dimension (WxHxD) (in / mm)	11.4 x 14.8 x 4.3 290 x 375 x 110	15.4 x 19.5 x 4.5 390 x 495 x 125	19.3 x 23.8 x 5.4 490 x 605 x 137
Net Weight (lb / Kg)	44.1 / 20	66.2 / 30	88.2 / 40
Electrical Requirement	100V - 240V 50/60Hz		
Electrical Plug Type	NEMA 5-15		
Power Consumption	60W		
Catalog Number	OS1511U	OS1512U	OS1513U

Incubated Orbital Shaker


Model	BIS-1	BIS-2	BIS-3
Maximum Load (lb / Kg)	4.4 / 2	8.8 / 4	19.8 / 9
Platform Size (W x D) (in / mm)	10.0 x 10.0 250 x 250	13.8 x 13.8 350 x 350	17.7 x 17.7 450 x 450
Speed Range (rpm)	40 - 250	40 - 300	
Frequency Accuracy (rpm)	±1		
Amplitude (mm)	20		
Control	LCD Microprocessor PID		
Heating Type	Forced-air Convection		
Temperature Range	Ambient + 9°F - 149°F / Ambient + 5°C - 65°C		
Display Resolution	0.1		
Temperature Uniformity (@37°C)	<±1.0		
Temperature Stability (@37°C)	±0.1		
Timer (hh:mm)	00:01 – 99:59		
Chamber Height (in / mm)	7.7 / 195	10.4 / 265	13.0 / 330
External Dimension (WxHxD) (in / mm)	15.4 x 23.2 x 14.6 390 x 590 x 370	19.3 x 27.2 x 17.7 490 x 690 x 450	23.2 x 32.5 x 21.7 590 x 825 x 550
Net Weight (lb / Kg)	70.6 / 32	92.6 / 42	114.6 / 52
Universal Clamp	Standard		
Electrical Requirement	120V/60Hz		
Electrical Plug Type	NEMA 5-15		
Power Consumption	450W	650W	1000W
Catalog Number	BS1511U	BS1512U	BS1513U

Shaker Accessories


Universal Platform

			Orbital Shaker			Incubated Shaker		
Accessory	SKU#	Dimension	BS-1	BS-2	BS-3	BIS-1	BIS-2	BIS-3
	P30016	250mm × 250mm	•			•		
	P30017	350mm × 350mm		•			•	
	P30018	450mm × 450mm			•			•

Universal Attachment


			Orbital Shaker			Incubated Shaker		
Accessory	SKU#	Dimension	BS-1	BS-2	BS-3	BIS-1	BIS-2	BIS-3
	P30004	250mm × 250mm x 85mm	•			•		
	P30005	350mm × 350mm x 100mm		•			•	
	P30006	450mm × 450mm x 100mm			•			•

Spring Wire Rack

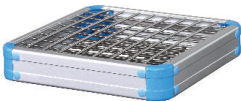
			Orbital Shaker			Incubated Shaker		
Accessory	SKU#	Dimension	BS-1	BS-2	BS-3	BIS-1	BIS-2	BIS-3
	P30001	250mm × 250mm	•			•		
	P30002	350mm × 350mm		•			•	
	P30003	450mm × 450mm			•			•

Mounting hardware included.

Cell Culture Spring Rack

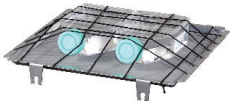
			Orbital Shaker			Incubated Shaker		
Accessory	SKU#	Dimension	BS-1	BS-2	BS-3	BIS-1	BIS-2	BIS-3
	P30007	250mm × 250mm	•			•		
	P30008	350mm × 350mm		•			•	
	P30009	450mm × 450mm			•			•

Rubber Mat

			Orbital Shaker			Incubated Shaker		
Accessory	SKU#	Dimension	BS-1	BS-2	BS-3	BIS-1	BIS-2	BIS-3
	P30013	250mm x 250mm	•			•		
	P30014	350mm x 350mm		•			•	
	P30015	450mm x 450mm			•			•


Mounting hardware included.

Flask Clamp

			Orbital Shaker			Incubated Shaker		
Accessory	SKU#	Flask Size	BS-1	BS-2	BS-3	BIS-1	BIS-2	BIS-3
	P33072	50ml	16	36	49	16	36	49
	P33073	100ml	9	23	36	9	23	36
	P33074	250ml	5	13	18	8	13	18
	P33075	500ml	4	8	16	4	8	16
	P33076	1000ml	–	4	9	–	4	9
	P33077	2000ml	–	1	4	–	–	4


Numbers shown indicate the total number of same-size flask clamps to fit on a universal platform. Flask clamps are sold individually. Requires purchase of a universal platform if needed. Visit our website to obtain the flask layout diagram.

Tube Rack

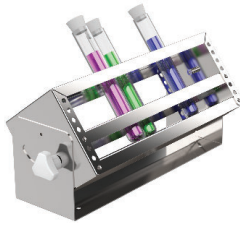		Tube Diameter x Number of Tubes per Rack	Orbital Shaker			Incubated Shaker		
Accessory	SKU#		BS-1	BS-2	BS-3	BIS-1	BIS-2	BIS-3
	P33078	Ø8mm x 105	2	3	4	2	3	4
	P33079	Ø10mm x 72						
	P33080	Ø12mm x 72						
	P33081	Ø14mm x 42						
	P33082	Ø17mm x 26						
	P33083	Ø25mm x 14						
	P33084	Ø30mm x 12						
	P33085	Ø35mm x 10						
	P33086	Ø50mm x 4						

Numbers shown indicate the total number of tube racks that fit on a universal platform. Tube racks are sold individually. Requires purchase of a universal platform if needed. Tube racks can be orientated from +45° to -45° and any angle in between.

Sticky Pad

			Orbital Shaker			Incubated Shaker		
Accessory	SKU#	Dimension	BS-1	BS-2	BS-3	BIS-1	BIS-2	BIS-3
	P30013	1450mm x 140mm	1	4	9	1	4	9

BEING's portfolio of laboratory equipment includes.

Incubators

BIF Series - Mechanical Convection Incubator
BIT Series - Natural Convection Incubator
BIC Series - Cooling Incubator

Ovens

BOF Series - Forced-air Drying Oven
BON Series - Natural Convection Drying Oven
BOV Series - Vacuum Oven

Shakers

BS Series - Orbital Shaker
BIS Series - Incubated Shaker

Stirrers

BMS Series - Square Plate Heated Magnetic Stirrer

Water Bath

BWB Series - General Purpose Water Bath
BWB Series - Dual Chamber Water Bath
BWS Series - Shaking Water Bath
BPC Series - Heat/Cooling Circulating Bath
BRC Series - Recirculating Chiller

Pumps

V Series - Diaphragm Pumps

Learn more at www.beinglab-usa.com
The 'Smart Choice' for laboratory equipment.

Your local dealer is.

novatech
INTERNATIONAL

sales@novatech-usa.com
www.novatech-usa.com
Tel: (866) 433-6682 Fax: (866) 433-6684
Tel: (281) 359-8538 Fax: (281) 359-0084

Shakers Brochure 03-2023 V1

Contact us

being

BEING Scientific Inc.

800 N. Haven Ave., Suite 428
Ontario, CA 91764

T: 800.278.1390
E: sales@beinglab-usa.com
www.beinglab-usa.com

Connect with us



©2023. All rights reserved.