

# Traceable® Fingertip Timer

**Easy-to-use, Place anywhere—  
Impressive Memory feature**

**Quick-to-program, three-key, countdown timer  
delivers error-free operation for any lab test**

- **Features:** count down, time in/out, alarm
- **Impressive memory remembers most-recently-programmed setting for recurring lab experiments**
- **Handy size** (1<sup>13</sup>/<sub>16</sub> x 3<sup>1</sup>/<sub>2</sub> x 1<sup>1</sup>/<sub>2</sub> inches)

## Features:

- Use anywhere: incubators, refrigerators, ovens, chambers, cabinets and walls
- Carry anywhere else with ease
- Individually serialized, calibrated, and certified Traceable® to NIST


### Traceable to NIST for accuracy

An individually-numbered Traceable® Certificate is provided which assures accuracy from an ISO/IEC 17025:2005 (1750.01) calibration laboratory accredited by A2LA. It indicates traceability to standards provided by NIST (National Institute of Standards and Technology).

## Alarm features:

- Repeat memory saves last programmed time, thus eliminating having to re-enter test times
- Once countdown is complete, an alarm will sound
- Manually turn alarm off or have the alarm automatically turn off after 30 seconds

Cat. No. 5042 Traceable® Fingertip Timer

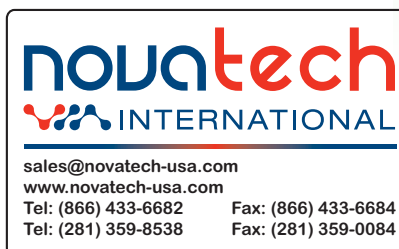

## Product Specifications

- ▶ **Timing Capacity:** 99 minutes, 59 seconds
- ▶ **Accuracy:** 0.01%
- ▶ **Resolution:** 1 second
- ▶ **Dimensions:** 1<sup>13</sup>/<sub>16</sub> x 3<sup>1</sup>/<sub>2</sub> x 1<sup>1</sup>/<sub>2</sub>"
- ▶ **Weight:** 1¼ ounces
- ▶ **Supplied:** battery, magnetic back, Velcro®, Traceable® Certificate