



B4000-16

- + Fully automatic; just press start
- + Rapid, complete sterilization
- + Extremely compact
- + Built in water tank, does not require connection to an external water supply

Despite the space-saving exterior dimensions of the BioClaves, the stainless steel sterilization chambers are surprisingly large, accommodating a variety of liquids, media, instruments, glassware, plasticware and other common research laboratory items. Both of these benchtop models include a mechanical and electrical safety interlock that prevents the door from being opened until the pressure has been released.




B4000-M

With large digital displays and fully automatic operation, all segments of the sterilization cycle (fill, sterilize, exhaust and dry) commence and run to completion by simply pressing the start button. An optional thermal mini-printer is also available for providing a record of sterilization parameters. In addition, a built in USB port allows for transfer of these parameters to a computer for printing or storage of cycle data.

Both BioClaves offer three cycle options (standard programs), which are preset for the sterilization of liquids, wrapped instruments or unwrapped instruments/plasticware. BioClaves are not designed for non-research applications such as medical diagnostics or foods. A fourth option "Dry Only" is also available for additional drying time to be added to the end of a cycle. For special applications, requiring variations in cycle parameters, all settings on the BioClaves are fully adjustable.



BioClave Mini

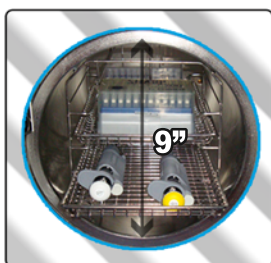

BioClave 16

Technical Data:

B4000-M (BioClave Mini)

Chamber Volume:	8L
Max. Pressure:	29 PSI / 2 Bar
Chamber Dimensions:	6.7 x 12.5 in. / 17 x 31cm
Tray Dimensions:	4.75 x 9.75 in. / 12 x 24.8 cm
Exterior Dimensions:	20 x 13.5 x 12.9 in.
<small>(W x D x H)</small>	51 x 34 x 33 cm
Weight:	74 lbs (33 kg)
Electrical:	115 or 230V, 50-60Hz, 950W
Warranty:	1 Year

B4000-16 (BioClave)

Chamber Volume:	16L
Max. Pressure:	29 PSI / 2 Bar
Chamber Dimensions:	9 x 13.75 in. / 23 x 35 cm
Tray Dimensions:	6.5 x 10.75 in. / 16.5 x 27 cm
Exterior Dimensions:	22 x 17.5 x 15.7 in.
<small>(W x D x H)</small>	56 x 44.5 x 40 cm
Weight:	89 lbs (40 kg)
Electrical:	115 or 230V, 50-60Hz, 1600W
Warranty:	1 Year

Ordering Information:

Item No.	Description:
B4000-M*	BioClave Mini Autoclave (8L)
B4000-16*	BioClave Autoclave (16L)
B4000-P*	Printer for BioClave series
B4000-PA	Extra roll of paper for B4000-P

*115V with US plug. For 230V (EU plug), please add (-E). This product is intended for research use only.