

GHS Compliant Wash Bottles



Multi-lingual Wide Mouth GHS Compliant Wash Bottles

- Multi-lingual wide mouth wash bottles featuring the patented DripLok® vapor venting valve that prevents solvent drips
- Indelibly printed in 4 languages (English, French, German and Spanish)
- Choose from four solvent labels, all including the new Globally Harmonized Symbols, NFPA diamond, risk phrases and safety information
- Fitted with color coded one-piece PP closures that maintain their shape in use
- Closures feature precision jet-tip aperture for fine control

LDPE MULTI-LINGUAL PRINTED SOLVENT VENTING WASH BOTTLES - 250ML

506485	-0001	-0002	-0003	-0004	-0005*
Capacity ml	250	250	250	250	250
Height mm	145	145	145	145	145
Diameter mm	60	60	60	60	60
Printed	Acetone	Ethanol	Isopropanol	Methanol	Set all 250ml
Pack Quantity	5	5	5	5	1

LDPE MULTI-LINGUAL PRINTED SOLVENT VENTING WASH BOTTLES - 500ML

506495	-0001	-0002	-0003	-0004	-0005*
Capacity ml	500	500	500	500	500
Height mm	165	165	165	165	165
Diameter mm	75	75	75	75	75
Printed	Acetone	Ethanol	Isopropanol	Methanol	Set all 500ml
Pack Quantity	5	5	5	5	1

Indelibly printed in 4 languages:
English, French, German & Spanish



*Sets include 1 of each printed bottle as well as the corresponding capacity custom write-on bottle.

Azlon Wash Bottles

Tracker Bottles, LDPE

- Allows full traceability of solutions
- All information is clearly visible to users, at all times
- Easy to re-use, time and time again
- DripLok® vapor venting version ensures safe use with volatile solvents
- Self-adhesive COSHH labels supplied with all bottles

With and without the Driplok® Vent these new LDPE wash bottles are designed to allow for full traceability of solutions. The bottle has been designed with a series of designated labeled areas, including spaces to record the following information.

- Originator
- Date of receipt
- Lot number
- Bottle contents
- First use date
- Expiration date
- Spaces for COSHH hazard symbols

506995	-0001	-0002	-0003	-0004	WGWPEN01
Description	Standard	Standard	DripLok®	DripLok®	Indelible
Capacity ml	250	500	250	500	Marker
Height mm	145	165	145	1165	Pen
Diameter mm	60	75	60	75	NA
Case Quantity	5	5	5	5	2




Custom Write-on Venting Wash Bottles



- Wide mouth for easy filling
- Eliminates unwanted jets of solvent
- Dispenses solvents efficiently
- Eliminates solvent dripping

These LDPE Driplok Wash Bottles with PP closures feature the only venting system that dispenses solvent as efficiently as a standard wash bottle. This venting system consists of a valve specifically designed to allow volatile solvent vapors to escape, but seals completely when in use to dispense a powerful, spatter-free jet of solvent. These easy to fill bottles are available blank with write on for customization.

506475	-0001	-0002	The Azlon® Driplok system (illustrated below) is the ONLY solvent vapor venting bottle that dispenses as efficiently as a normal wash bottle due to the patented ball valve system. No spit...no sputter!
Capacity ml	250	500	
Height mm	145	165	
Diameter mm	60	75	
Printed	Write-on	Write-on	
Pack Quantity	5	5	

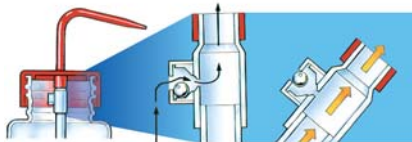

NEW Wash Bottle, Autoclavable, PP



- Autoclavable at 121°C
- PP 500mL bottle fitted with one-piece PP closure (34mm neck ID)
- Closures feature precision jet-tip aperture for fine control
- Dimensions 163 height x 75mm diameter

506985



Biological Stain Bottle Kit

Stain series Driplok venting wash bottles preprinted for five popular biological stains. Printed for Lugol's Iodine Solution, Crystal Violet, Gram Color De-coloring Solution, Gram Iodine Solution and Safranin Solution. Printed information includes hazard warning symbols where relevant and instructions for use. Dark colored stain bottles are printed in light colors for easy recognition of contents. Consists of five assorted 250ml bottles labeled as above.

507005 Biological Stain Set



175 Humboldt Street, Suite 300 Rochester, NY 14610
Telephone 800-334-7585 | Fax 585-334-0241

www.dynalon.com

Dynalon®
Labware
quality . value . service