

- **VersaCap®** technology provides the flexibility of interchangeable connectors
- **EZclear®** body features advanced ergonomic design
- **EZfil™** wide-mouth design for easy filling and cleaning
- **VersaBarb®** technology provides interchangeable spigot options
- **Made in the USA** 

Caution

- Do not use carboys under pressure or full vacuum
- Do not use carboys at temperatures greater than the rated maximum
- Do not store oxidizing agents in carboys
- Do not flame carboys
- Do not use carboy spigot as a handle
- Do not expose to organic compounds. Carboy will become prone to stress fracturing


Foxx Life Sciences' **EZclear®** Carboy is an advanced fluid handling system developed with the end-user in mind. The ergonomically designed **EZclear®** body is coupled with our **VersaCap®** technology to provide the user an easy-to-use and easy-to-clean carboy.

The Copolyester Carboy is made from a clear, BPA-free, USP Class VI, FDA Grade resin. This Carboy is suitable for aqueous solutions, as well as most life sciences buffer solutions.


The **EZclear®** Carboy's body includes design features that allows the user to easily hold and manipulate during lifting, filling, and pouring. The **EZclear®** Carboy features include large top handles, molded hand grips in the base, easy-to-read metric graduation marks certified to $\pm 5\%$ accuracy, and molded-in material identification markings.

The Copolyester 20L Carboy offers an **EZfil™** 120mm neck to simplify filling and cleaning, while all other sizes feature an 80mm neck. All carboys are compatible with our **VersaCap®** technology for interchangeable top connections.

ISO 9001
ISO 13485

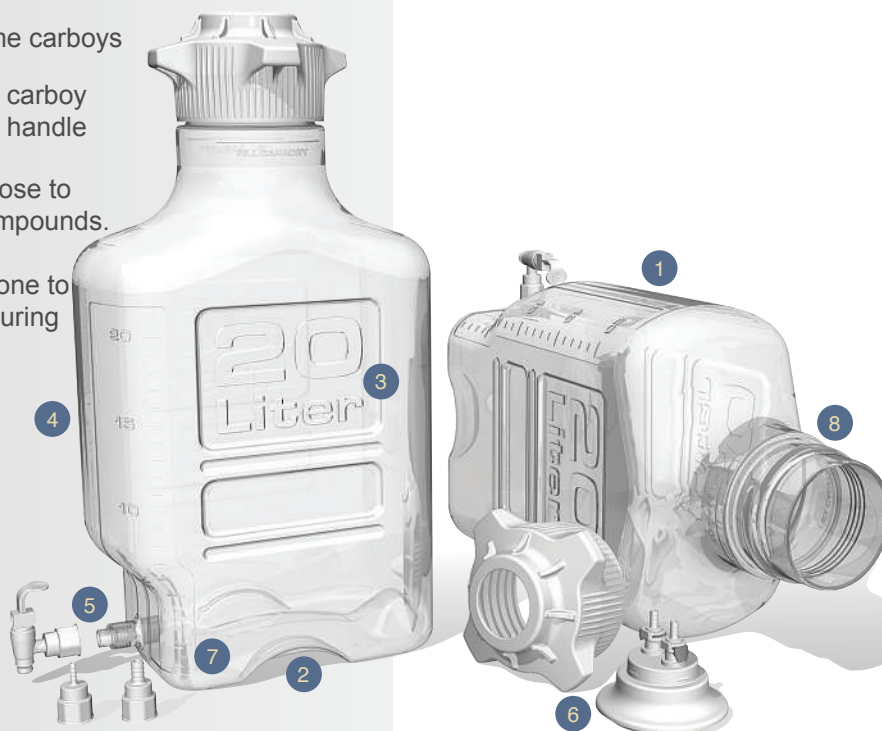

- 1 Novel shape maximizes bench and shelf space
- 2 Bottom handles provide a secure grip while lifting and transporting
- 3 Molded-in carboy identification marks
- 4 Large, clear graduation marks with $\pm 5\%$ accuracy
- 5 Autoclavable spigot and **VersaBarb®** technology
- 6 **VersaCap®** technology with versatile adapters and replaceable gaskets
- 7 **EZclear®** body design allows for maximum solution recovery
- 8 **EZfil™** wide-mouth shape simplifies filling and cleaning

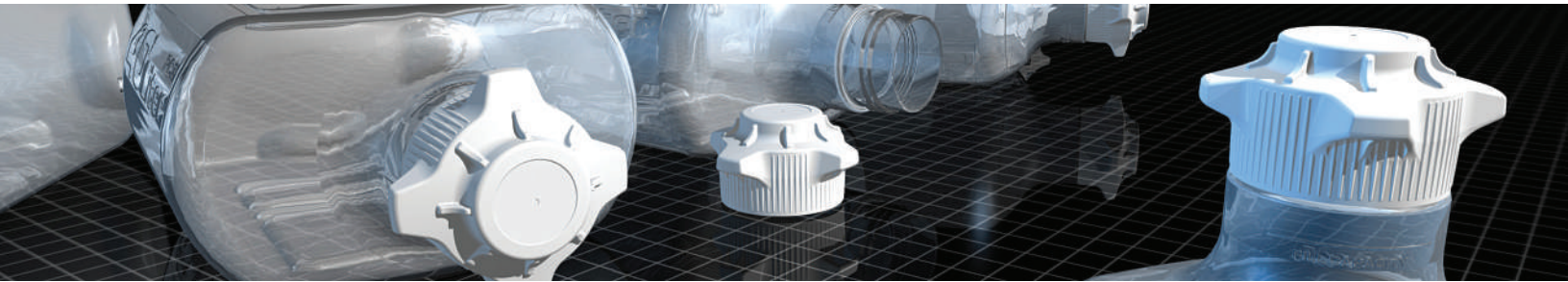

Product Number	Description	Capacity (L)	Material	Spigot	Cap Size (mm)	Qty/Pack
157-0111-OEM	Carboy, CP, 2.5L, VersaCap® 80	2.5	CP	No	80	1
157-1111-OEM	Carboy, CP, 5L, VersaCap® 80	5	CP	No	80	1
157-2111-OEM	Carboy, CP, 10L, VersaCap® 80	10	CP	No	80	1
157-3121-OEM	Carboy, CP, 20L, VersaCap® 120	20	CP	No	120	1



Wide range of caps and adapters to fit all your application needs.
See **VersaCap®**, **VersaBarb®** and Spigot Technology data sheet.

novatech
INTERNATIONAL

sales@novatech-usa.com
www.novatech-usa.com
Tel: (866) 433-6682 Fax: (866) 433-6684
Tel: (281) 359-8538 Fax: (281) 359-0084

We recommend testing for your specific application before use.