

- **VersaCap®** technology provides the flexibility of interchangeable connectors
- **EZgrip™** body features advanced ergonomic design
- **EZfil™** features a wide-mouth design for easy filling and cleaning
- **VersaBarb®** technology provides interchangeable spigot options

- **Made in the USA**



### Caution

- Do not use carboys under pressure or full vacuum
- Do not use carboys at temperatures greater than the rated maximum
- Do not store oxidizing agents in carboys
- Do not flame carboy
- Do not use carboy spigot as a handle


Foxx Life Sciences' Polypropylene **EZgrip™** Carboy is an advanced fluid handling system developed with the end-user in mind. The ergonomically designed **EZgrip™** body is coupled with **VersaCap®** technology to provide the user with an easy-to-use and easy-to-clean carboy.

The **EZgrip™** Polypropylene Carboy is made from translucent USP Class VI, FDA Grade autoclavable resin, which offers good chemical resistance and heat stability. The recommended autoclaving cycle is 15 PSI for 20 minutes at 121°C.


The Carboy body includes design features that allow the user to easily hold and manipulate during lifting, filling, and pouring, while **VersaCap®** technology provides interchangeable, leak-proof top connections. The **EZgrip™** Carboy features include large top handles, molded hand grips in the base, easy-to-read metric graduation marks certified to  $\pm 5\%$  accuracy, and the carboy material molded into the side for easy identification.

The 20L, 40L, and 75L carboys offer an **EZfil™** 120mm neck to simplify filling and cleaning, while all other sizes feature an 83mm neck. All Carboys are compatible with our **VersaCap®** technology for interchangeable top connections.

ISO 13485

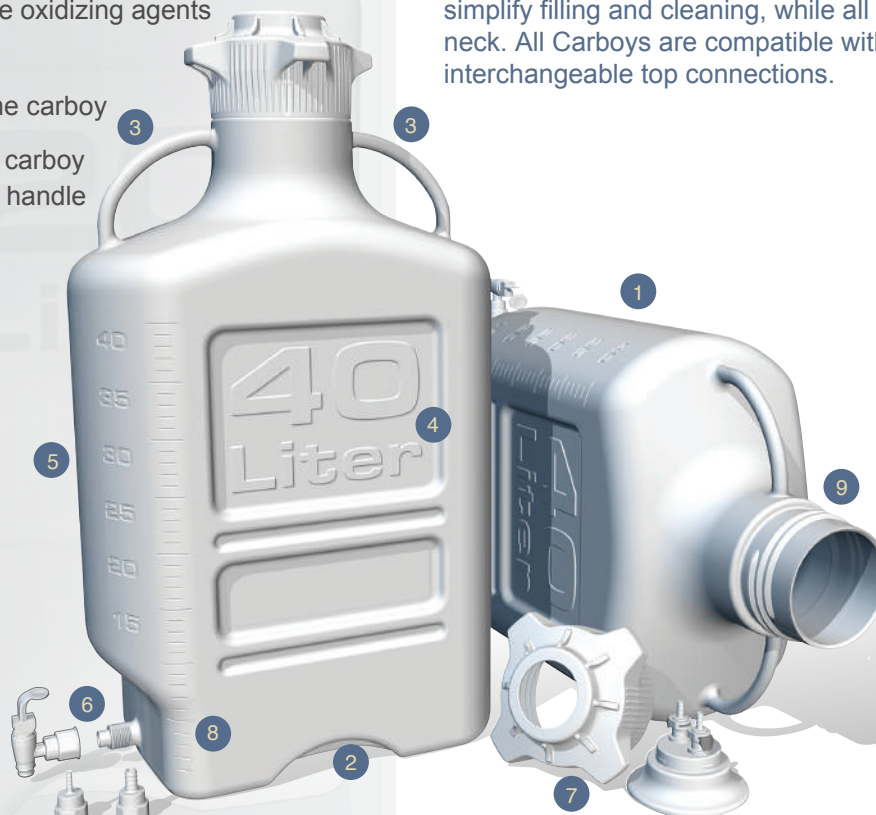

- 1 Novel shape maximizes bench and shelf space
- 2 Bottom handles provide a secure grip while lifting and transporting
- 3 Large top handles for a full hand grip
- 4 Molded carboy identification marks
- 5 Large, clear graduation marks with  $\pm 5\%$  accuracy
- 6 Autoclavable spigot and **VersaBarb®** technology
- 7 **VersaCap®** technology with versatile adapters and replaceable gasket
- 8 **EZgrip™** body design allows for maximum solution recovery
- 9 **EZgrip™** wide-mouth shape for easy filling and cleaning

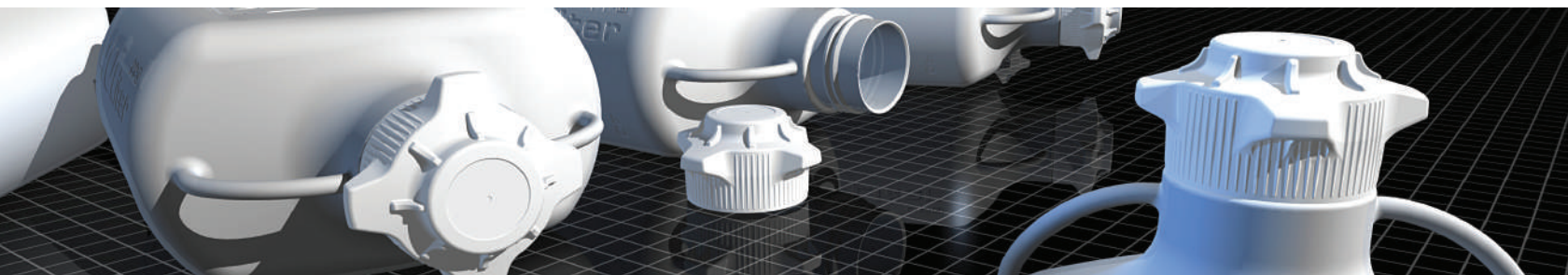

Product Number	Description	Capacity (L)	Material	Spigot	Cap Size (mm)	Qty/Pack
155-0111-OEM	Carboy, PP, 2.5L, <b>VersaCap®</b> 83B	2.5	PP	No	83	1
155-1111-OEM	Carboy, PP, 5L, <b>VersaCap®</b> 83B	5	PP	No	83	1
155-1211-OEM	Carboy, PP, 5L, Spigot, <b>VersaCap®</b> 83B	5	PP	Yes	83	1
155-2111-OEM	Carboy, PP, 10L, <b>VersaCap®</b> 83B	10	PP	No	83	1
155-2211-OEM	Carboy, PP, 10L, Spigot, <b>VersaCap®</b> 83B	10	PP	Yes	83	1
155-2131-OEM	Carboy, PP, 10L, 3" Sanitary Neck	10	PP	No	N/A	1
155-3111-OEM	Carboy, PP, 20L, <b>VersaCap®</b> 83B	20	PP	No	83	1
155-3211-OEM	Carboy, PP, 20L, Spigot, <b>VersaCap®</b> 83B	20	PP	Yes	83	1
155-3131-OEM	Carboy, PP, 20L, 3" Sanitary Neck	20	PP	No	N/A	1
155-3121-OEM	Carboy, PP, 20L, <b>VersaCap®</b> 120mm	20	PP	No	120	1
155-3221-OEM	Carboy, PP, 20L, Spigot, <b>VersaCap®</b> 120mm	20	PP	Yes	120	1
155-5121-OEM	Carboy, PP, 40L, <b>VersaCap®</b> 120mm	40	PP	No	120	1
155-5221-OEM	Carboy, PP, 40L, Spigot, <b>VersaCap®</b> 120mm	40	PP	Yes	120	1
155-5131-OEM	Carboy, PP, 40L, 3" Sanitary Neck	40	PP	No	N/A	1
155-7121-OEM	Carboy, PP, 75L, No Spigot, <b>VersaCap®</b> 120mm	75	PP	No	120	1



BPA-Free



USP Class VI

**novatech**  
INTERNATIONAL

sales@novatech-usa.com

www.novatech-usa.com

Tel: (866) 433-6682

Fax: (866) 433-6684

Tel: (281) 359-8538

Fax: (281) 359-0084