

LABSTRONG REPLACEMENT FOR BARNSTEAD 740883 AND 740895




The LabStrong 740883LS replacement heating element is comparable to the Thermo Scientific™ Barnstead™ 740883 and 740895 replacement heating elements. The LabStrong 740883LS is an economical and reliable replacement heating element for Thermo Scientific Barnstead Mega-Pure™ Glass Stills. The 240V 2500w heating element is suitable for operation in the MP-3A (Model numbers A440367 and A440696), MP-6A (Model numbers A440518 and A440697), MP-11A (Model numbers A440118 and A440117) and MP-12A (Model numbers A442011 and A442012) Thermo Scientific Barnstead Mega Pure Glass Stills.

The LabStrong heating element offers a new and improved design that is directly compatible to the Thermo Scientific Barnstead part numbers 740883 and 740895. The heating element converts electricity into heat through the process of Joule heating. Electric current through the element encounters resistance, resulting in heating of the element. The heating of the element causes the water in the boiler of the glass still to boil. This heating allows for the process of distillation to occur. The heating element provides outstanding quality at a lower cost than the competitor's.

Product Specifications:

- Part Number: 740883LS
- Comparable Thermo Scientific Barnstead Part Number: 740883 and 740895
- Product Type: Replacement Heating Element for Water Purification Systems
- Compatible With: Thermo Scientific Barnstead Mega-Pure Water Stills (MP-3A, MP-6A, MP-11A an MP-12A)
- Electrical Voltage: 240 V
- Power: 2500 watts
- Current: 10.41 amps per element
- Resistance: 21.9 to 24.3 OHMS
- Compliance: RoHS Compliant

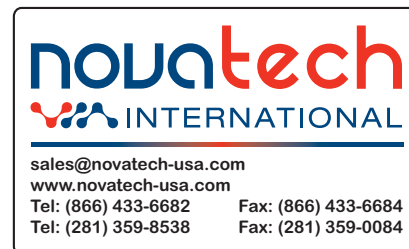

- Made in USA with Global Materials -