

Tygon® 2475 Tubing

Plasticizer Free Thermoplastic Tubing

Due to the absence of plasticizers, the Tygon® 2475 thermoplastic tubing can withstand aggressive acids, bases, and solvents that would react with or extract plasticizer from a standard PVC tube. The clear tubing allows you to visually inspect the fluid path to ensure consistency of flow. Designed with lower gas permeation than standard silicone tubing, sensitive fluids are protected from moisture loss and oxidation.

The surface smoothness of the tubing helps reduce potential microbial buildup, critical in bioprocess applications. The tubing also features a very low absorption of aqueous substances, minimizing the risk of fluid alteration in single or repeat use applications.

Materials

Thermoplastic Elastomer (TPE)

- Good chemical resistance

TPE compounds based on SBS and/or SEBS exhibit good resistance to water and various solutions. However, some degradation in the compounds can be seen at long-term exposure to organic solvents, oils, and fuels.



CERTIFICATION

Cole-Parmer will provide a certificate of compliance FREE at your request where applicable.

APPLICATIONS

- Aseptic filling and dispensing systems, diagnostic equipment, nuclear equipment, lab analytical instrumentation, cosmetic production, and cell and tissue culture transport

FEATURES/BENEFITS

- Low gas permeation
- Smooth inner surface inhibits particle entrapment
- Free from added plasticizers or oils
- Meets USP Class VI criteria
- Ideal for general-purpose, laboratory, and single-use applications
- Free technical application support available to assist with tubing selection


USA	+1.800.323.4340
	+1.847.549.7600
Canada	+1.800.363.5900
China	86.21.5109.9909
France	+33 (0) 1486 37800

Germany	+49 (0) 9377 92030
India	+1.800.266.1244
Italy	+39.02.84349215
UK	+44 (0) 1480.272279
All others	+1.847.549.7600

Specification and Ordering Table

Dimensions: in. (mm)			Max pressure: PSI (bar)	Item number	Reel length (ft)
Tubing ID	Tubing OD	Wall			
1/16"	3/16"	1/16" (1.6)	85 (5.8)	50113-72	50
1/8"	1/4"	1/16" (1.6)	50 (3.4)	50113-73	50
3/16"	5/16"	1/16" (1.6)	40 (2.7)	50113-74	50
1/4"	3/8"	1/16" (1.6)	30 (2.06)	50113-75	50
5/16"	7/16"	1/16" (1.6)	18 (1.2)	50113-76	50
3/8"	5/8"	1/16" (1.6)	20 (1.4)	50113-77	50
1/2"	3/4"	1/8" (3.2)	29 (2)	50113-78	50
5/8"	7/8"	1/8" (3.2)	25 (1.7)	50113-79	50
3/4"	1-1/4"	1/8" (3.2)	21 (1.4)	50113-80	50
1"	1-3/8"	3/16" (4.8)	20 (1.4)	50113-81	50

TYGON® 2475 TUBING

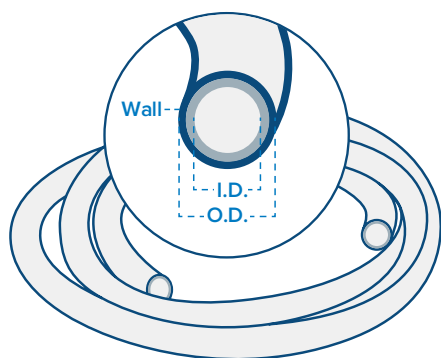

novatech
INTERNATIONAL

sales@novatech-usa.com

www.novatech-usa.com

Tel: (866) 433-6682

Fax: (866) 433-6684

Tel: (281) 359-8538

Fax: (281) 359-0084

Cole-Parmer®
an antylia scientific company